


LCP Anterolateral Distal Tibia Plate 3.5

Product No.	Hole	Direction
22574109	9H	Right
22574111	11H	Right
22574113	13H	Right
22574115	15H	Right
22574117	17H	Right
22574119	19H	Right
22574209	9H	Left
22574211	11H	Left
22574213	13H	Left
22574215	15H	Left
22574217	17H	Left
22574219	19H	Left

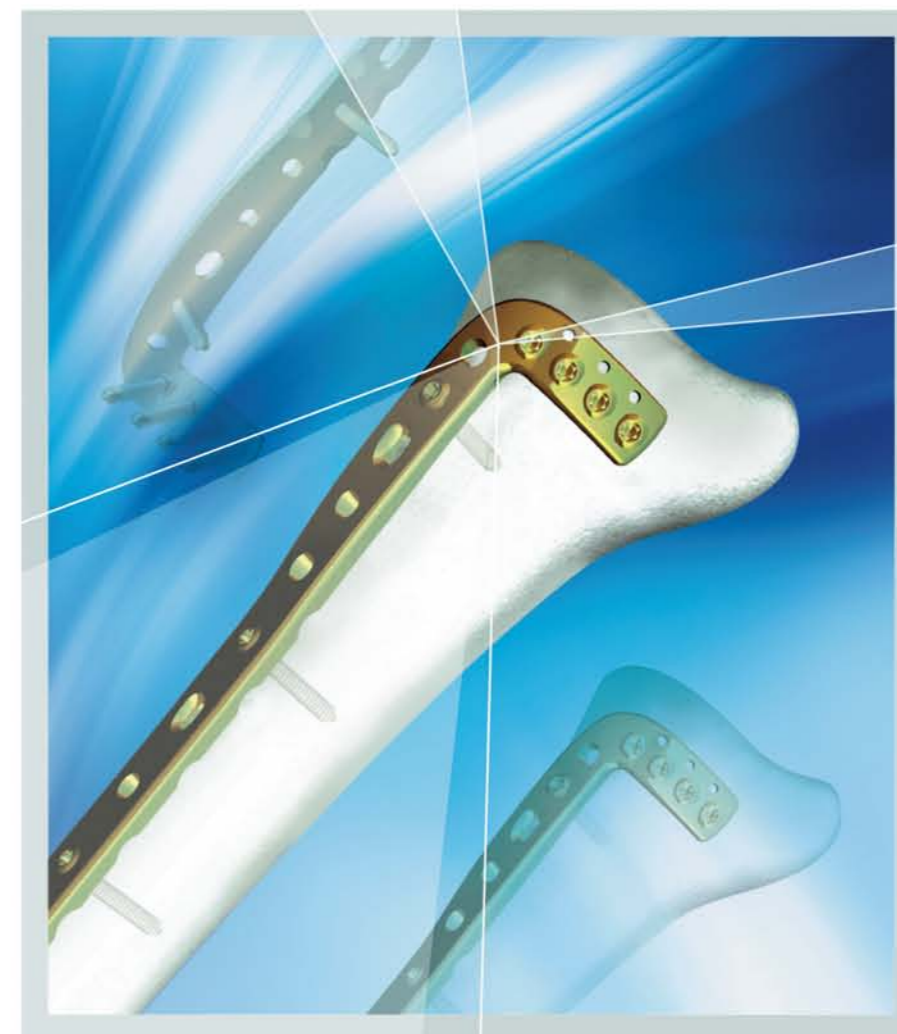

Monoloc

LCP Anterolateral Distal Tibia L Plate 3.5

CHANGZHOU KANGHUI MEDICAL INNOVATION CO., LTD

Add: No.11 North Changjiang Road, Xinbei Zone, Changzhou, Jiangsu 213022, P.R.China
Tel: +86-519-85139851 Fax: +86-519-85128628
www.kanghui.com
Shanghai Office
Add: Building 18, No. 1000, Jinhai Road, PuDong District, Shanghai, China 201206
Tel: +86-21-50319916 Fax: +86-21-50312913

T - 201407

Key Features & Benefits

- Precontoured to fit the anterolateral distal tibia.
- Particularly indicated for osteopenic bone.
- Tapered tip for sub muscular insertion.
- Distal locking screws provide support for the articular surface.
- Left and right plates.



Indications

- The 3.5mm LCP Anterolateral Distal Tibia Plate is indicated for fractures, osteotomies and nonunions of the distal tibia, particularly in osteopenic bone.

Design features

- A. Three K - wire holes in the head, parallel to the joint, accept 2.0 mm K - wires to temporarily fix the plate to the distal tibia, temporarily reduce articular fragments, and show proximity to the joint.
- B. Four distal head holes angle 7 ° inferiorly to capture the posterior malleolus.
- C. 3 Elongated holes aids in plate positioning and compression.
- D. Tapered tip for easy insertion.

